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Dear customer,

We are very pleased that you have decided in favor of a product from our range. Our name stands for high quality and thoroughly tested products in the fields of heat, gentle therapy, blood pressure/diagnosis, weight, massage and air. Please read these instructions for use carefully and comply with the directions.

With best regards
Your BEURER Team

Items included in the package

- Callus remover
- 2x Coarse roller
- 1x Fine roller
- Protective cap for nail dust
- Cleaning brush
- 2x 1.5 V AA batteries
- These instructions for use

1. Important notes

Read through these instructions for use carefully before first use. Keep them to hand for later use and ensure that they are accessible to other users.

IMPORTANT:

Please ensure that no hair can get caught in the rotating attachments when using the device. Use a hair tie for your safety.

During use, a considerable amount of heat may be generated on the skin.

- Check results regularly during treatment. This is particularly advisable in the case of diabetics, as they are less sensitive to pain and are liable to become injured more easily. Always use with extreme care. If in doubt, ask your doctor.
- This instrument must be used only for the purpose for which it has been developed, and in the manner described in the instructions for use. Any improper use can be dangerous.
- In case of prolonged intensive use (e.g. milling of dead skin from your feet) the device might get very hot. In order to prevent skin burns in these situations, there should be longer breaks between each application. For your own safety, constantly monitor heat build up in the device. This particularly applies to people who are insensitive to heat.
- The device is only for private use and is not intended for medical or commercial purposes.
- For reasons of hygiene, the device should only be used on one person.
- The manager is not responsible for damage or injury caused by improper or incorrect use.

- Ensure before use that there is no visible damage to the instrument and accessories. In the event of doubt, do not use them and contact your dealer or the customer service address stated.
- Keep children away from the packaging materials. Risk of suffocation.
- The appliance must be checked frequently for signs of wear or damage. If such signs are present, or if the appliance has been used improperly, it must be taken to either the manufacturer or the dealer before further use.
- Switch off the instrument immediately in the event of defects or malfunctions.
- Repairs may be carried out only by the customer service or authorized dealers.
- Before making a claim, please check the batteries first and replace them if necessary.
- On no account try to repair the instrument yourself!
- Use the instrument only with the accessories provided.
- Never operate the appliance without supervision, particularly when children are in the vicinity.
- Do not use on small animals/animals.
- Prevent all contact with water (except when cleaning with a slightly moistened cloth!). Water must never enter the inside of the appliance. Never immerse the appliance in water. Never use the device in the bath, in the shower, in a swimming pool, above a sink filled with water or in an area with extremely high humidity (e.g. the bathroom).
- Keep the instrument away from sources of heat.
- Do not use the instrument under blankets, pillows etc.
- Do not use the device for periods longer than 15 minutes without a pause. After this time, wait an interval of at least 15 minutes to prevent the motor from overheating.
- Hands must always be dry when using the device.
- This device may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the consequent risks of use.
- Do not use the device on areas of sensitive or irritated skin.
- Children should be supervised to ensure that they do not use the device as a toy.
- Cleaning and user maintenance must not be performed by children unless supervised.

Notes on handling batteries

- If your skin or eyes come into contact with battery fluid, flush out the affected areas with water and seek medical assistance.
-  Choking hazard! Small children may swallow and choke on batteries. Store the batteries out of the reach of small children.
- Observe the plus (+) and minus (-) polarity signs.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- Protect the batteries from excessive heat.
-  Risk of explosion! Never throw batteries into a fire.
- Do not charge or short-circuit batteries.

- If the device is not to be used for a long period, take the batteries out of the battery compartment.
- Use identical or equivalent battery types only.
- Always replace all batteries at the same time.
- Do not use rechargeable batteries.
- Do not disassemble, split or crush the batteries.

2. Preparation

In order to use the Callus remover, you first need to insert the supplied batteries. To insert the batteries, follow the steps below:

1. Remove the battery compartment lid on the rear of the device.
2. Insert the two supplied 1.5 V AA (alkaline type LR6) batteries. Make sure that the batteries are inserted with the correct polarity, according to the label. Do not use rechargeable batteries.
3. Close the battery compartment lid carefully.

3. Application


3.1 General

- This device is intended for the treatment of feet only (pedicure).
- To change an attachment, press the unlock button and take the attachment off the device.
- If you wish to remove highly callused areas, use the coarse callus attachment.
- If you wish to remove slightly callused areas or to smooth your skin, use the fine callus attachment.
- Switch on the device by moving the sliding switch forwards. When doing so, also press the switch-on guard on the sliding switch downwards.
- The sandpaper attachments are characterised by their extremely high effectiveness.
- Do not exert heavy pressure and always guide the attachments carefully over the surface to be treated.
- Move the instrument under light pressure in circular movements slowly over the parts to be treated.
- Note that the cutting attachments will have a reduced effect if the skin is softened or moist. For this reason, do not carry out any pretreatment in a bath.
- Do not remove all the callus skin in order to retain the natural protection of the skin.
- After every use, apply moisturizing cream to the parts treated.

Changing the roller:

Proceed as follows to change the roller:

- Hold the callus remover in your hand as shown.
- With your other hand, push the roller to the left and pull it up and off.

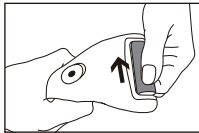4. Cleaning/care

- Clean the instrument with a slightly moistened cloth. You can also moisten the cloth with mild soap suds to remove more stubborn dirt.
- The sandpaper attachments must only be cleaned using a dry cloth or brush.
- Do not drop the instrument.
- Make sure that no water enters the instrument. But if this should happen, do not use the instrument again until it is completely dry.
- Do not use chemical detergent or scouring agents for cleaning.
- Never immerse the device in water or other liquids.

5. Replacing the battery

If the inserted batteries are empty, you must replace them.

To replace the batteries, follow the instructions in the "Initial use" chapter.

6. Replacement parts and wearing parts

Replacement part	Item number and/or order number
Sandpaper attachment replacement set	573.07

Commercially available or from the corresponding service address (according to the service address list)

7. Disposal

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life. Dispose of the unit at a suitable local collection or recycling point. If you have any questions, please contact the local authorities responsible for waste disposal.

The empty, completely flat batteries must be disposed of through specially designated collection boxes, recycling points or electronics retailers. You are legally required to dispose of the batteries.

The codes below are printed on batteries containing harmful substances:

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